Data sheet

dNTP Set

Cat. No N0098

4X40 µmol (4X0.4 mL of 100 mM solution)

Description

dNTP set includes 4 vials, a separate vial of each dNTP (dATP, dCTP, dGTP, dTTP). Each dNTP are in aqueous solution at 100 mM.

Applications

For use in all molecular biology applications, including PCR, real-time PCR, high fidelity and long PCR, cDNA synthesis, RT-PCR, and DNA sequencing.


The synthesis of a complementary DNA/cDNA strand by enzymatic polymerization of deoxynucleotide triphosphates (dNTPs) is the basic principle of DNA amplification, sequencing and cDNA synthesis techniques.

TruePurity dNTPs Certification:

- ✓ Functionally tested in 18 kb long range PCR (template dilution from 20 to 1000 pg).
- Not contamination with bacterial and human DNA.
- ✓ No detectable activity of DNAse, protease or phosphatase.

Storage: Store at -20°C, short term (up to one week) exposure to ambient temperature is possible.

dATP

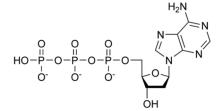

Formula: $C_{10}H_{16}N_5O_{12}P_3$ Molecular weight: 491.18 g·mol-1 Concentration: 100 mM

dCTP

Formula: $C_9H_{16}N_3O_{13}P_3$

Molecular weight: 467.15 g·mol-1 Concentration: 100 mM

dGTP

Formula: C₁₀H₁₆N₅O₁₃P₃ Molecular weight: 507.18 g·mol-1

Concentration: 100 mM

dTTP

Formula: C₁₀H₁₇N₂O₁₄P₃

Molecular weight: 482.17 g⋅mol-1

Concentration: 100 mM

PRODUCT USE LIMITATION

This product is developed, designed and sold exclusively for research purposes and in vitro use only. The product was not tested for use in diagnostics or for drug development, nor is it suitable for administration to humans or animals. Please refer to www.canvaxbiotech.com for Material Safety Data Sheet of the product.



Tallaght Business Park Whitestown, Dublin 24, Ireland D24 RFK3

Tel: (01) 4523432 Fax: (01) 4523967 Web: www.labunlimited.com Quatro House, Frimley Road, Camberley, United Kingdom GU16 7ER

Tel: 08452 30 40 30 Fax: 08452 30 50 30 E-mail: info@labunlimited.com E-mail: info@labunlimited.co.uk Web: www.labunlimited.co.uk





Notes



Tallaght Business Park Whitestown, Dublin 24, Ireland D24 RFK3

Tel: (01) 4523432 Web: www.labunlimited.com

Quatro House, Frimley Road, Camberley, United Kingdom GU16 7ER

Fax: (01) 4523967 Fax: 08452 30 50 30
E-mail: info@labunlimited.com E-mail: info@labunlimited.co.uk Web: www.labunlimited.co.uk

